


simpass CTO

PTCA BALLOON CATHETER

A Reinvented **simpass** for **CTO** interventions

Renewed Distal Shaft
Tapered Down Design &
Hybrid Polymer Material

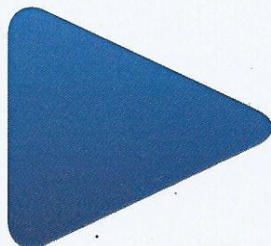

0,016" Entry Profile
With Atraumatic Soft Tip

Low Profile Balloon that
folded with SMART Folding
Technique

Specially Designed Tapered Down and Hybrid Polymer Distal Shaft

Tapered-Down:


In Distal Part, 2,5 F Profile For Deliverability
In Proximal Part, 2,7 F Pusher Body Segment For Superior Pushability

Hybrid Polymer:

In Distal Part, Soft Polymer Material For Superlative Flexibility
In Proximal Part, Stiff Polymer Material for Extra Support

Low Entry Profile by SMART Folding Technique

0,016" Low Profiled Distal Tip For First Contact With Lesion

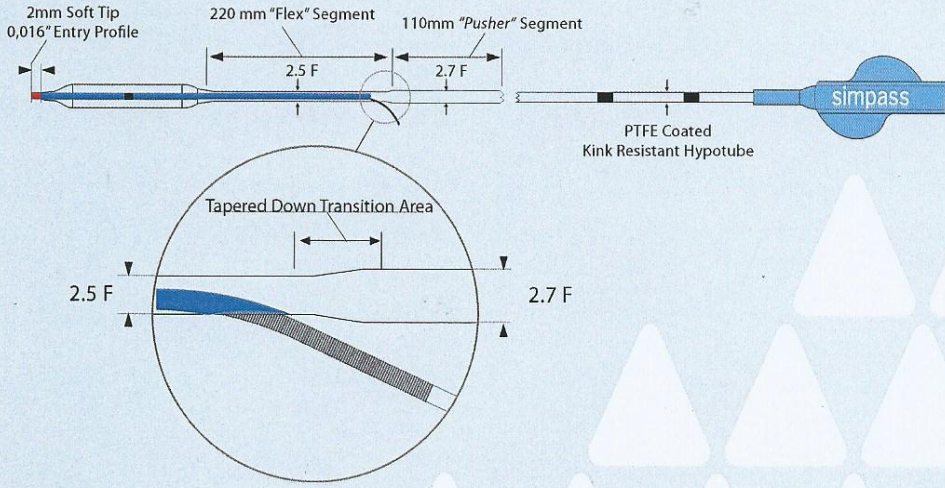

- ▶ Conical shaped tip for CTO lesions
- ▶ 0.016\" Entry profile
- ▶ Smooth lesion crossing with single marker
- ▶ 2.5 F Flexible distal shaft
- ▶ Pusher segment between distal and proximal shaft for better pushability
- ▶ Proximal shaft with increased kink resistance
- ▶ Wide range of product sizes

Balloon Diameters (mm)	Balloon Lengths (mm)		
	10	15	20
1.00	ST010010	ST010015	ST010020
1.25	ST012510	ST012515	ST012520
1.50	ST015010	ST015015	ST015020
2.00	ST020010	ST020015	ST020020
2.50	ST025010	ST025015	ST025020

Pressure (atm)	Balloon Diameters (mm)				
	1.00	1.25	1.50	2.00	2.50
4	0.90	1.17	1.41	1.90	2.39
6	0.95	1.22	1.47	1.97	2.47
8*	1.00	1.25	1.50	2.01	2.53
10	1.03	1.28	1.53	2.05	2.56
12	1.06	1.31	1.56	2.08	2.59
14	1.09	1.34	1.58	2.11	2.63
16**	1.12	1.37	1.60	2.13	2.66
18	1.17	1.41	1.63	2.16	2.69
20	1.20	1.45	1.66	2.20	2.73

simeks

SIMEKS Tıbbi Ürünler Sanayi ve Ticaret Limited Şirketi
 İstanbul Endüstri ve Ticaret Serbest Bölgesi Orjin Cad.
 10. Sok. No: 113/3 34953 Tuzla / İSTANBUL
 Tel: +90(216)3940071 | Faks: +90(216)3940081
 Web: <http://www.simeksmedical.com>

MADE IN TURKEY / TÜRK MALIDIR